


Category - Endodontic

FLEXER®

ONE

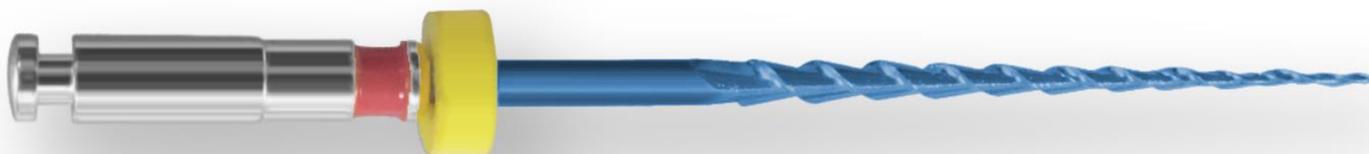

Features

- Equipped with Twin Alloys [IonTech 6.0 Niti + Monel Alloy U400].
- Equipped with Twin Technology [S.A.W & Heat Activation].
- Non-Cutting Tip - Avoids Perforation.
- Variable Cross Section.

Description

FLEXER - ONE, unique design makes it easier to bend the endodontic files to follow the anatomy of the canal during treatments, without deforming it and without forming sharp edges. The product allows the professionals to perform the procedures precisely and safely.

Key specifications

- Single File system.
- Long durability.
- High Aggressive Cutting Efficiency.
- Continuous clockwise Rotation Function.
- 300 – 500% Improved Fracture Resistance.
- Accurate shaping to protect the Dentin effectively.

Packaging

- Pack of 3 pieces.

SIZE	T1
TIP	#25
TAPER (mm)	0.06

Application Sequence

- * Initially use hand file till #25.
- * Use Mx (or) Gates drill for orifice opening.
- * Use glide path.
- * Use **FO** till working length, to shape and complete the canals.